



CONFORT EXCEPTIONAL^{1,2} ȘI STABILITATE^{1,3} IDEALĂ CU LENTILE TORICE DIN GAMA TOTAL

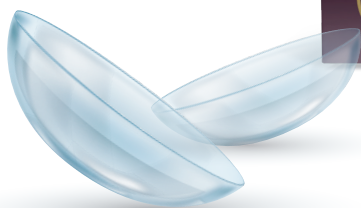

2
CARACTERISTICI
COMUNE

&
2
TEHNOLOGII
UNICE


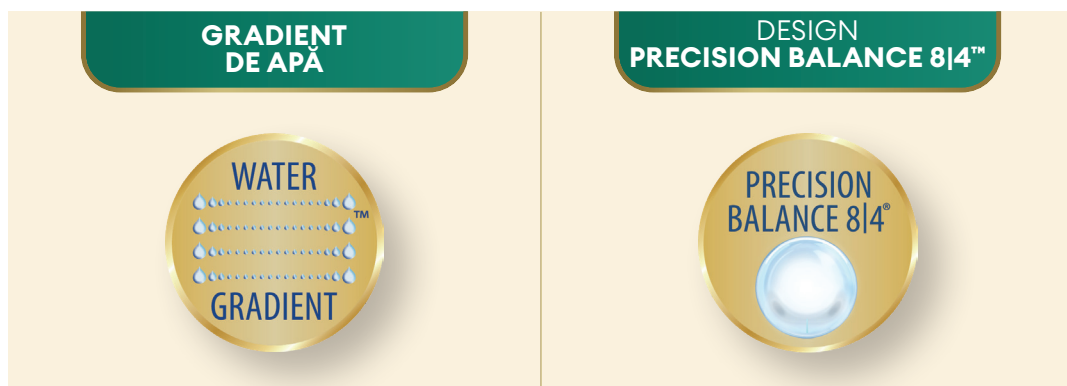


Referințe: 1. In a clinical trial to evaluate on-eye performance of TOTAL30® for Astigmatism lenses where n=69; Alcon data on file, 2021. 2. In a clinical trial to assess overall performance of DAILIES TOTAL1® for Astigmatism lenses where n=134 patients; Alcon data on file, 2021. 3. In a clinical trial to evaluate stability of axis orientation of DAILIES TOTAL1® for Astigmatism lenses where n=47; Alcon data on file, 2020.

Alcon

2 CARACTERISTICI COMUNE

GRADIENT DE APĂ

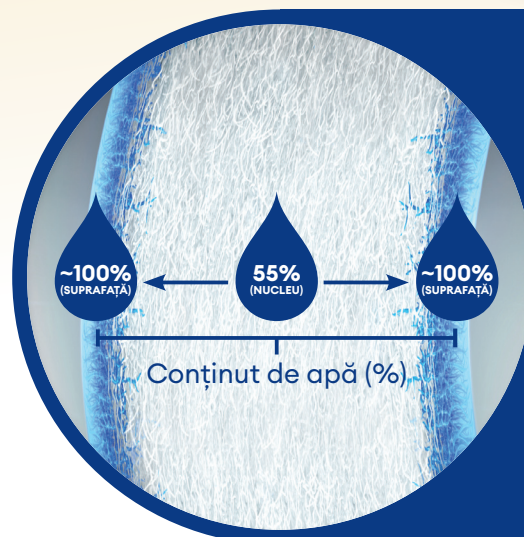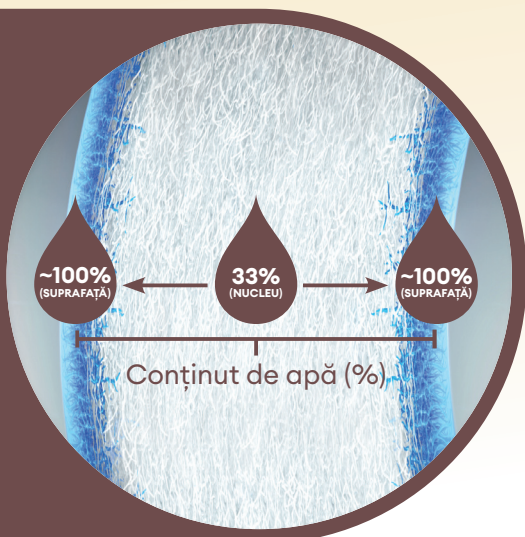

Ambele lentile oferă o combinație ideală de transmisibilitate a oxigenului și conținut ridicat de apă la suprafața exterioară a lentilei, **astfel încât tot ceea ce atinge ochiul este o pernă moale de umiditate.**

PRECISION BALANCE 8|4™ DESIGN



STABILITATE EXCELENTĂ^{5,6}

Oferă o vedere mai bună de 20/20 cu o stabilitate excelentă.*⁷

* În cazul pacienților care nu au ambliopie.

CONFORT IDEAL^{6,7}

Lentila este mai subțire în zona marchajului de la ora 6 minimizând astfel interacțiunea cu pleopa.^{5,6}


Un singur
marcaj la
ora 6

≤5°
oscilație la
clipire^{5,6}

100%
DINTRE UTILIZATORI
AU EXPERIMENTAT
O OSCILAȚIE LA
CLIPIRE DE ≤5°^{5,6}

STABILIZARE ÎN
~60
SECUENDE^{5,6}

≥ **95%**
RATĂ DE SUCCES A
PRIMULUI FITTING^{5,6*}

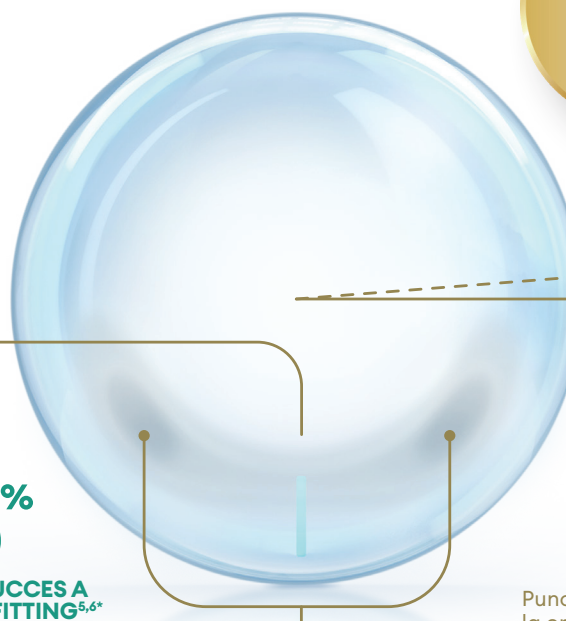


Puncte de stabilizare
la ora 8 și ora 4 care
reduce interacțiunea
cu pleopa inferioară

*Based on lens movement, centration, and rotation at initial fitting.


Referințe: 1. In vitro analysis of leihfilcon A contact lenses outermost surface softness and correlation with water content; Alcon data on file, 2021. 2. In vitro analysis of lens oxygen permeability, water content, and surface imaging; Alcon data on file, 2021. 3. Angelini TE, Nixon RM, Dunn AC, et al. Viscoelasticity and mesh-size at the surface of hydrogels characterized with microrheology. Invest Ophthalmol Vis Sci. 2013;54:E-abstract 500. 4. Thekveli S, Qui Y, Kapoor Y, et al. Structure-property relationship of deilefilcon A lenses. Cont Lens Anterior Eye. 2012;35(Suppl 1):e14. 5. In a clinical trial to evaluate stability of axis orientation of DAILIES TOTAL1® for Astigmatism lenses where n=47; Alcon data on file, 2020. 6. In a clinical trial to evaluate on-eye performance of TOTAL30® for Astigmatism lenses where n=69; Alcon data on file, 2021. 7. In a clinical trial to assess overall performance of DAILIES TOTAL1® for Astigmatism lenses where n=134 patients; Alcon data on file, 2021.

TRANSMISIBILITATE EXCELENTĂ A OXIGENULUI¹⁻³


DAILIES TOTAL1® ASTIGMATISM

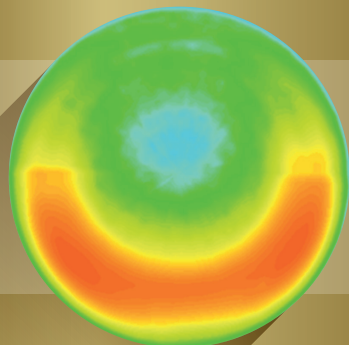


**Cel mai mare Dk/t pe întreaga suprafață a
lentilei, inclusiv zonele periferice^{2*}**

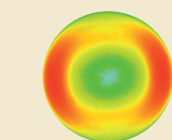

ACUVUE® OASYS 1-DAY
for Astigmatism
with hydraluxe

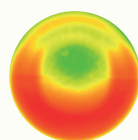

MyDay® Toric

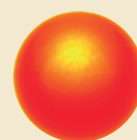


clariti® 1 day Toric



1-DAY ACUVUE® MOIST
for Astigmatism
with Lacreon





Biotrue® ONEday
for Astigmatism

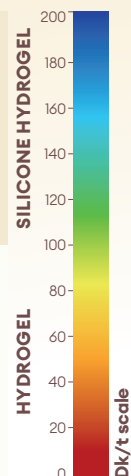


**DAILIES TOTAL1®
ASTIGMATISM
ESTE SUPERIOR
LENTILELOR CONCURENTE²**


TOTAL30® ASTIGMATISM

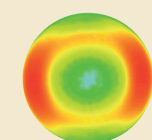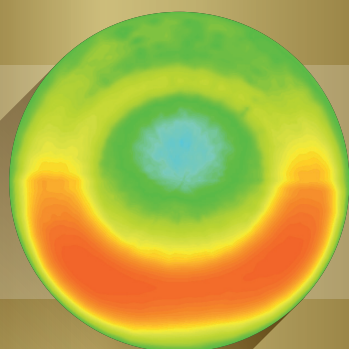


ACUVUE® OASYS
for Astigmatism with
HydraClear Plus®

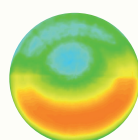

Biofinity® Toric

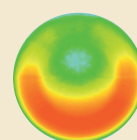

Bausch + Lomb ULTRA®
for Astigmatism

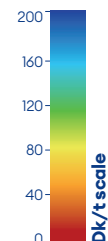

**TOTAL30® ASTIGMATISM ARE
CEA MAI MARE RATĂ DE
RESPIRABILITATE COMPARTIV CU
ALTE LENTILE TORICE LUNARE^{1,3**}**

*Dk/t = 127 (@-3.00/-1.25 x180)


**Dk/t = 112 (@-3.00/-1.25 x 180)

*Trademarks are the property of their respective owners.


Referințe: 1. Laboratory measurement of oxygen permeability by polarographic method; Alcon data on file, 2019. 2. Based on in vitro measurement of thickness profiles on unworn lenses; Alcon data on file, 2019. 3. Based on in vitro measurements of oxygen transmissibility and oxygen permeability values of competitive reusable toric contact lenses; Alcon data on file, 2018.

TIMP SUPERIOR DE RUPERE A STRATULUI DE UMIDITATE^{1-5*}

Datorită gradientului de apă, DAILIES TOTAL1® și TOTAL30® au
TIMP ÎNDELUNGAT DE RUPERE A UMIDITAȚII DE LA SUPRAFAȚA LENTILEI^{1-5*}

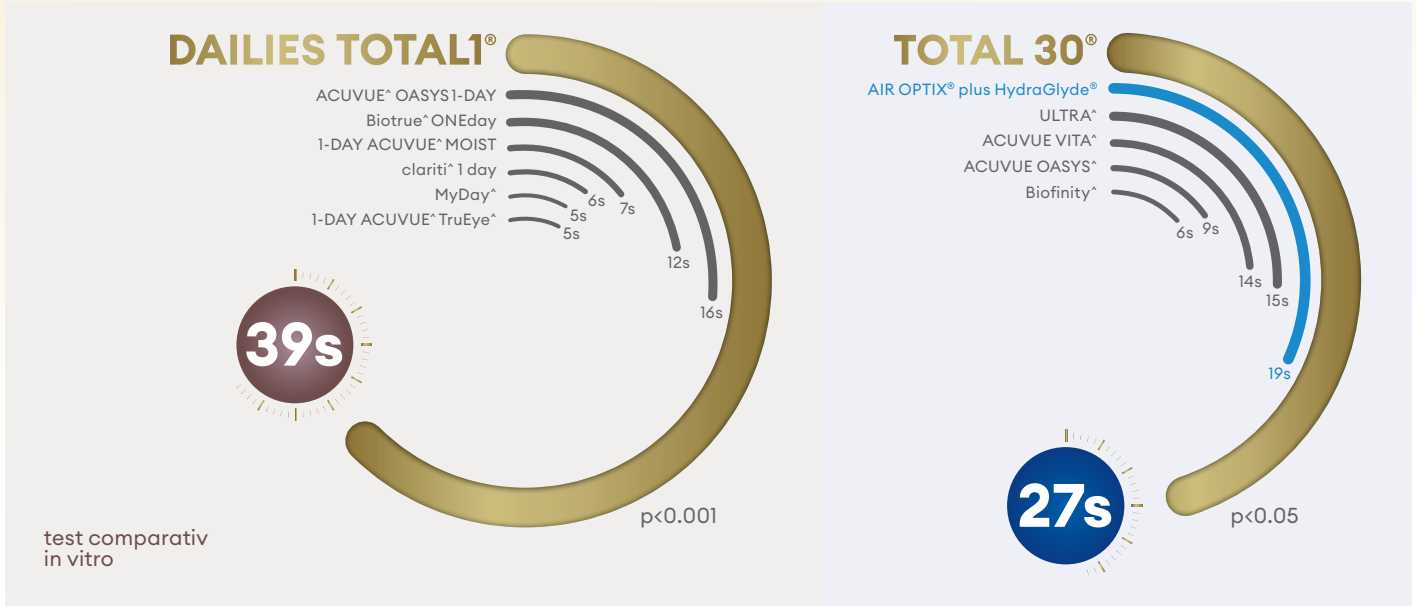

INFORMAȚII PRODUS

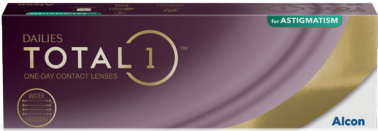

Material	delefilcon A			lehfilcon A		
Rază de curbură (mm)	8.6			8.6		
Grosime centrală (@-3.00D/-1.25 x 180, mm)	0.11			0.10		
Modulus (MPa)	0.7			0.6		
Diametru (mm)	14.5			14.5		
Parametrii	Sferă	Cilindru	Ax	Sferă	Cilindru	Ax
	+4.00D to +0.25D (0.25D pași)	-0.75D, -1.25D, -1.75D	10°, 20°, 70°, 80°, 90°, 100°, 110°, 160°, 170°, 180°	+6.00 to -6.00 (0.25D pași)	-0.75 to -2.25 (0.50D pași)	10° to 180° (10° pași)
		-2.25D	10°, 20°, 160°, 170°, 180°			
	Plano to -6.00D (0.25D pași)	-0.75D, -1.25D, -1.75D	în sensul acelor de ceasornic în pași de 10°	+6.00 to +8.00 și -6.50 to -10.00 (0.50D pași)		
		-2.25D	10°, 20°, 70°, 80°, 90°, 100°, 110°, 160°, 170°, 180°			
	-6.50D to -8.00D (0.50D pași)	-0.75D, -1.25D, -1.75D	10°, 20°, 70°, 80°, 90°, 100°, 110°, 160°, 170°, 180°			
		-2.25D	10°, 20°, 160°, 170°, 180°			
Conținut de apă la suprafață	Aproximativ 100%			Aproximativ 100%		
Conținut de apă în centru	33%			55%		
Dk/t	127 @-3.00D/-1.25 x 180			123 @-3.00D/-1.25 x 180		
Ambalaj	5buc. (lentile de probă), 30buc., 90 buc.			1 buc. (probă), 3 buc., 6 buc		
Filtre	filtru UV clasa I, filtru HEVL†			filtru UV clasa I, filtru HEVL†		

*Based on an in vitro study of sphere lenses wherein wettability was measured using the iDDrop System. All lenses were tested in an identical manner, soaked in a PBS (phosphatebuffered saline solution) for 16 hours +/- 2 hours.
**UV absorbing contact lenses are NOT substitutes for protective UV absorbing eyewear, such as UV absorbing goggles or sunglasses because they do not completely cover the eye and surrounding area. The patient should continue to use UV absorbing eyewear as directed.
† There is no demonstrated clinical benefit to a 34% reduction in visible light at wavelengths below 450 nm.
^Trademarks are the property of their respective owners.
Referințe: 1. Based on in vitro studies of unworn TOTAL30® and competitive sphere lenses (Biofinity, Ultra, ACUVUE OASYS and ACUVUE VITA) using average iDDrop water break-up time (WBUT); Alcon data on file, 2021. 2. Based on an in vitro study on unworn DAILIES TOTAL1® and competitive sphere lenses (1-DAY ACUVUE MOIST, MyDay and clariti 1day) using average iDDrop water break-up time (WBUT); Alcon data on file, 2019. 3. Based on an in vitro study on unworn ACUVUE OASYS 1-DAY and ACUVUE TrueEye sphere lenses using average iDDrop water break-up time (WBUT); Alcon data on file, 2019. 4. Based on an in vitro study on unworn Infuse sphere lenses testing average iDDrop water break-up time (WBUT) and surface softness; Alcon data on file, 2020. 5. Based on an in vitro study on unworn BioTrue ONEday sphere lenses using average iDDrop water break-up time (WBUT); Alcon data on file, 2018.

GHID DE ADAPTARE

INSTRUCȚIUNI PENTRU UTILIZAREA TABELULUI DE CONVERSIE

NOTĂ: Nu este necesar să faceți dumneavoastră conversia distanței vertex deoarece tabelul are valorile deja calculate!

- 1

Identificați puterea sferică pentru rețeta de ochelari în partea stângă, sub denumirea "dioptrie sferică ochelari"
- 2

Identificați lentilele de contact de probă DAILIES TOTAL1® sau TOTAL30® pentru Astigmatism cu semnul (-) la puterea cilindrului de la ochelarii purtătorului sau va trebui să faceți transpoziție.
- 3

Pentru lentile de contact de probă, selectați DAILIES TOTAL1® sau TOTAL30® pentru Astigmatism apoi alegeți axul cu cea mai apropiată valoare de valoarea refracției sau de axul prescris pentru ochelari.
- 4

Atunci când comandați lentilele de contact finale, trebuie să țineți cont de rotația lentilelor de contact de probă ca fiind corespunzătoare.



ACOPERĂ

94%
DIN ASTIGMAȚI

Corecție vertex: 12mm

DIOPTRIE CILINDRICĂ OCHELARI (CYL)

DIOPTRIE SFERICĂ OCHELARI (SPH)

PARAMETRII LENTILĂ DE CONTACT DE PROBĂ

	-0.75		-1.00		-1.25		-1.50		-1.75		-2.00		-2.25		-2.50		-2.75		-3.00	
	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl
-9.00	-8.00	sph	-8.00	-0.75	-8.00	-0.75	-8.00	-1.25	-8.00	-1.25	-8.00	-1.75	-8.00	-1.75	-8.00	-1.75	-8.00	-2.25	-8.00	-2.25
-8.75	-8.00	sph	-8.00	-0.75	-8.00	-0.75	-8.00	-1.25	-8.00	-1.25	-8.00	-1.75	-8.00	-1.75	-8.00	-1.75	-8.00	-2.25	-8.00	-2.25
-8.50	-8.00	sph	-7.50	-0.75	-7.50	-0.75	-7.50	-1.25	-7.50	-1.25	-7.50	-1.75	-7.50	-1.75	-7.50	-1.75	-7.50	-2.25	-7.50	-2.25
-8.25	-7.75	sph	-7.50	-0.75	-7.50	-0.75	-7.50	-1.25	-7.50	-1.25	-7.50	-1.75	-7.50	-1.75	-7.50	-1.75	-7.50	-2.25	-7.50	-2.25
-8.00	-7.75	sph	-7.50	-0.75	-7.50	-0.75	-7.50	-1.25	-7.50	-1.25	-7.50	-1.75	-7.50	-1.75	-7.50	-1.75	-7.50	-2.25	-7.50	-2.25
-7.75	-7.50	sph	-7.00	-0.75	-7.00	-0.75	-7.00	-1.25	-7.00	-1.25	-7.00	-1.75	-7.00	-1.75	-7.00	-1.75	-7.00	-2.25	-7.00	-2.25
-7.50	-7.25	sph	-7.00	-0.75	-7.00	-0.75	-7.00	-1.25	-7.00	-1.25	-7.00	-1.75	-7.00	-1.75	-7.00	-1.75	-7.00	-2.25	-7.00	-2.25
-7.25	-7.00	sph	-6.50	-0.75	-6.50	-0.75	-6.50	-1.25	-6.50	-1.25	-6.50	-1.75	-6.50	-1.75	-6.50	-1.75	-6.50	-2.25	-6.50	-2.25
-7.00	-6.75	sph	-6.50	-0.75	-6.50	-0.75	-6.50	-1.25	-6.50	-1.25	-6.50	-1.75	-6.50	-1.75	-6.50	-1.75	-6.50	-2.25	-6.50	-2.25
-6.75	-6.50	sph	-6.00	-0.75	-6.00	-0.75	-6.00	-1.25	-6.00	-1.25	-6.00	-1.75	-6.00	-1.75	-6.00	-1.75	-6.00	-2.25	-6.00	-2.25
-6.50	-6.00	-0.75	-6.00	-0.75	-6.00	-0.75	-6.00	-1.25	-6.00	-1.25	-6.00	-1.75	-6.00	-1.75	-6.00	-1.75	-6.00	-2.25	-6.00	-2.25
-6.25	-5.75	-0.75	-5.75	-0.75	-5.75	-0.75	-5.75	-1.25	-5.75	-1.25	-5.75	-1.75	-5.75	-1.75	-5.75	-1.75	-5.75	-2.25	-5.75	-2.25
-6.00	-5.50	-0.75	-5.50	-0.75	-5.50	-0.75	-5.50	-1.25	-5.50	-1.25	-5.50	-1.75	-5.50	-1.75	-5.50	-1.75	-5.50	-2.25	-5.50	-2.25
-5.75	-5.50	-0.75	-5.50	-0.75	-5.50	-0.75	-5.50	-1.25	-5.50	-1.25	-5.50	-1.75	-5.50	-1.75	-5.50	-1.75	-5.50	-2.25	-5.50	-2.25
-5.50	-5.25	-0.75	-5.25	-0.75	-5.25	-0.75	-5.25	-1.25	-5.25	-1.25	-5.25	-1.75	-5.25	-1.75	-5.25	-1.75	-5.25	-2.25	-5.00	-2.25
-5.25	-5.00	-0.75	-5.00	-0.75	-5.00	-0.75	-5.00	-1.25	-5.00	-1.25	-5.00	-1.75	-5.00	-1.75	-5.00	-1.75	-5.00	-2.25	-5.00	-2.25
-5.00	-4.75	-0.75	-4.75	-0.75	-4.75	-0.75	-4.75	-1.25	-4.75	-1.25	-4.75	-1.75	-4.75	-1.75	-5.00	-1.75	-4.75	-2.25	-4.75	-2.25
-4.75	-4.50	-0.75	-4.50	-0.75	-4.50	-0.75	-4.50	-1.25	-4.50	-1.25	-4.50	-1.75	-4.50	-1.75	-4.75	-1.75	-4.50	-2.25	-4.50	-2.25
-4.50	-4.25	-0.75	-4.25	-0.75	-4.50	-0.75	-4.25	-1.25	-4.25	-1.25	-4.25	-1.75	-4.25	-1.75	-4.50	-1.75	-4.25	-2.25	-4.25	-2.25
-4.25	-4.00	-0.75	-4.00	-0.75	-4.25	-0.75	-4.00	-1.25	-4.00	-1.25	-4.00	-1.75	-4.00	-1.75	-4.25	-1.75	-4.00	-2.25	-4.25	-2.25
-4.00	-4.00	-0.75	-4.00	-0.75	-4.00	-0.75	-4.00	-1.25	-4.00	-1.25	-4.00	-1.75	-4.00	-1.75	-4.25	-1.75	-4.00	-2.25	-4.00	-2.25
-3.75	-3.75	-0.75	-3.75	-0.75	-3.75	-0.75	-3.75	-1.25	-3.75	-1.25	-3.75	-1.75	-3.75	-1.75	-3.75	-1.75	-3.75	-2.25	-4.00	-2.25
-3.50	-3.50	-0.75	-3.50	-0.75	-3.50	-0.75	-3.50	-1.25	-3.50	-1.25	-3.50	-1.75	-3.50	-1.75	-3.50	-1.75	-3.50	-2.25	-3.75	-2.25
-3.25	-3.25	-0.75	-3.25	-0.75	-3.25	-0.75	-3.25	-1.25	-3.25	-1.25	-3.25	-1.75	-3.25	-1.75	-3.25	-1.75	-3.25	-2.25	-3.50	-2.25
-3.00	-3.00	-0.75	-3.00	-0.75	-3.00	-0.75	-3.00	-1.25	-3.00	-1.25	-3.00	-1.75	-3.00	-1.75	-3.00	-1.75	-3.00	-2.25	-3.25	-2.25
-2.75	-2.75	-0.75	-2.75	-0.75	-2.75	-1.25	-2.75	-1.25	-2.75	-1.25	-2.75	-1.75	-2.75	-1.75	-2.75	-1.75	-2.75	-2.25	-3.00	-2.25
-2.50	-2.50	-0.75	-2.50	-0.75	-2.50	-1.25	-2.50	-1.25	-2.50	-1.25	-2.50	-1.75	-2.50	-1.75	-2.75	-1.75	-2.50	-2.25	-2.75	-2.25
-2.25	-2.25	-0.75	-2.25	-0.75	-2.25	-1.25	-2.25	-1.25	-2.25	-1.75	-2.25	-1.75	-2.25	-1.75	-2.50	-1.75	-2.25	-2.25	-2.50	-2.25
-2.00	-2.00	-0.75	-2.00	-0.75	-2.00	-1.25	-2.00	-1.25	-2.00	-1.75	-2.00	-1.75	-2.00	-1.75	-2.25	-1.75	-2.00	-2.25	-2.25	-2.25
-1.75	-1.75	-0.75	-1.75	-0.75	-1.75	-1.25	-1.75	-1.25	-1.75	-1.75	-1.75	-1.75	-2.00	-1.75	-2.00	-1.75	-1.75	-2.25	-2.00	-2.25
-1.50	-1.50	-0.75	-1.50	-0.75	-1.50	-1.25	-1.50	-1.25	-1.50	-1.75	-1.50	-1.75	-1.75	-1.75	-1.50	-2.25	-1.50	-2.25	-1.75	-2.25
-1.25	-1.25	-0.75	-1.25	-0.75	-1.25	-1.25	-1.25	-1.25	-1.25	-1.75	-1.25	-1.75	-1.50	-1.75	-1.25	-2.25	-1.50	-2.25	-1.50	-2.25
-1.00	-1.00	-0.75	-1.00	-0.75	-1.00	-1.25	-1.00	-1.25	-1.00	-1.75	-1.00	-1.75	-1.25	-1.75	-1.00	-2.25	-1.25	-2.25	-1.25	-2.25
-0.75	-0.75	-0.75	-0.75	-0.75	-0.75	-1.25	-0.75	-1.25	-0.75	-1.75	-0.75	-1.75	-1.00	-1.75	-0.75	-2.25	-1.00	-2.25	-1.00	-2.25
-0.50	-0.50	-0.75	-0.50	-0.75	-0.50	-1.25	-0.50	-1.25	-0.50	-1.75	-0.50	-1.75	-0.75	-1.75	-0.50	-2.25	-0.75	-2.25	-0.75	-2.25
-0.25	-0.25	-0.75	-0.25	-0.75	-0.25	-1.25	-0.25	-1.25	-0.25	-1.75	-0.25	-1.75	-0.50	-1.75	-0.25	-2.25	-0.50	-2.25	-0.50	-2.25
+0.00	+0.00	-0.75	+0.00	-0.75	+0.00	-1.25	+0.00	-1.25	+0.00	-1.75	+0.00	-1.75	-0.25	-1.75	+0.00	-2.25	-0.25	-2.25	-0.25	-2.25
+0.25	+0.25	-0.75	+0.25	-0.75	+0.25	-1.25	+0.25	-1.25	+0.25	-1.75	+0.25	-1.75	0.00	-1.75	+0.25	-2.25	+0.00	-2.25	+0.00	-2.25
+0.50	+0.50	-0.75	+0.50	-0.75	+0.50	-1.25	+0.50	-1.25	+0.50	-1.75	+0.50	-1.75	+0.25	-1.75	+0.50	-2.25	+0.25	-2.25	+0.25	-2.25
+0.75	+0.75	-0.75	+0.75	-0.75	+0.75	-1.25	+0.75	-1.25	+0.75	-1.75	+0.75	-1.75	+0.50	-1.75	+0.75	-2.25	+0.50	-2.25	+0.50	-2.25
+1.00	+1.00	-0.75	+1.00	-0.75	+1.00	-1.25	+1.00	-1.25	+1.00	-1.75	+1.00	-1.75	+0.75	-1.75	+1.00	-2.25	+0.75	-2.25	+0.75	-2.25
+1.25	+1.25	-0.75	+1.25	-0.75	+1.25	-1.25	+1.25	-1.25	+1.25	-1.75	+1.25	-1.75	+1.00	-1.75	+1.25	-2.25	+1.00	-2.25	+1.00	-2.25
+1.50	+1.50	-0.75	+1.50	-0.75	+1.50	-1.25	+1.50	-1.25	+1.50	-1.75	+1.50	-1.75	+1.25	-1.75	+1.50	-2.25	+1.25	-2.25	+1.25	-2.25
+1.75	+1.75	-0.75	+1.75	-0.75	+1.75	-1.25	+1.75	-1.25	+1.75	-1.75	+1.75	-1.75	+1.50	-1.75	+1.75	-2.25	+1.50	-2.25	+1.50	-2.25
+2.00	+2.00	-0.75	+2.00	-0.75	+2.00	-1.25	+2.00	-1.25	+2.00	-1.75	+2.00	-1.75	+1.75	-1.75	+2.00	-2.25	+1.75	-2.25	+1.75	-2.25
+2.25	+2.25	-0.75	+2.25	-0.75	+2.25	-1.25	+2.25	-1.25	+2.25	-1.75	+2.25	-1.75	+1.75	-1.75	+2.25	-2.25	+2.00	-2.25	+2.00	-2.25
+2.50	+2.50	-0.75	+2.50	-0.75	+2.50	-1.25	+2.50	-1.25	+2.50	-1.75	+2.50	-1.75	+2.25	-1.75	+2.50	-2.25	+2.25	-2.25	+2.25	-2.25
+2.75	+2.75	-0.75	+2.75	-0.75	+2.75	-1.25	+2.75	-1.25	+2.75	-1.75	+2.75	-1.75	+2.50	-1.75	+2.75	-2.25	+2.50	-2.25	+2.50	-2.25
+3.00	+3.00	-0.75	+3.00	-0.75	+3.00	-1.25	+3.00	-1.25	+3.00	-1.75	+3.00	-1.75	+2.75	-1.75	+3.00	-2.25	+2.75	-2.25	+2.75	-2.25
+3.25	+3.25	-0.75	+3.25	-0.75	+3.25	-1.25	+3.25	-1.25	+3.25	-1.75	+3.25	-1.75	+3.00	-1.75	+3.25	-2.25	+3.00	-2.25	+2.75	-2.25
+3.50	+3.50	-0.75	+3.50	-0.75	+3.50	-1.25	+3.50	-1.25	+3.50	-1.75	+3.50	-1.75	+3.25	-1.75	+3.50	-2.25	+3.25	-2.25	+3.00	-2.25
+3.75	+3.75	-0.75	+3.75	-0.75	+3.75	-1.25	+3.75	-1.25	+3.75	-1.75	+3.75	-1.75	+3.50	-1.75	+3.75	-2.25	+3.50	-2.25	+3.25	-2.25
+4.00	+4.00	-0.75	+4.00	-0.75	+4.00	-1.25	+4.00	-1.25	+4.00	-1.75	+4.00	-1.75	+3.75	-1.75	+4.00	-2.25	+4.00	-2.25	+3.75	-2.25
+4.25																			+4.00	-2.25
	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cy	sph	cyl	sph	cyl	sph	cyl

PARAMETRII LENTILĂ DE CONTACT DE PROBĂ

* acoperire pentru pacienții astigmati cu cilindru de -2.50 sau mai mic.



DIOPTRIE CILINDRICĂ OCHELARI (CYL)

	-0.75		-1.00		-1.25		-1.50		-1.75		-2.00		-2.25		-2.50		-2.75		-3.00		-3.25		-3.50	
	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl	sph	cyl
-11.50	-10.00	sph	-10.00	-0.75	-10.00	-0.75	-10.00	-0.75	-10.00	-1.25	-10.00	-1.25	-10.00	-1.25	-10.00	-1.75	-10.00	-1.75	-10.00	-2.25	-10.00	-2.25	-10.00	-2.25
-11.25	-10.00	sph	-10.00	-0.75	-10.00	-0.75	-10.00	-0.75	-10.00	-1.25	-10.00	-1.25	-10.00	-1.25	-10.00	-1.75	-10.00	-1.75	-10.00	-2.25	-10.00	-2.25	-10.00	-2.25
-11.00	-10.00	sph	-9.50	-0.75	-9.50	-0.75	-9.50	-0.75	-9.50	-1.25	-9.50	-1.25	-9.50	-1.25	-9.50	-1.75	-9.50	-1.75	-9.50	-2.25	-9.50	-2.25	-9.50	-2.25
-10.75	-9.50	sph	-9.50	-0.75	-9.50	-0.75	-9.50	-0.75	-9.50	-1.25	-9.50	-1.25	-9.50	-1.25	-9.50	-1.75	-9.50	-1.75	-9.50	-2.25	-9.50	-2.25	-9.50	-2.25
-10.50	-9.50	sph	-9.00	-0.75	-9.00	-0.75	-9.00	-0.75	-9.00	-1.25	-9.00	-1.25	-9.00	-1.25	-9.00	-1.75	-9.00	-1.75	-9.00	-2.25	-9.00	-2.25	-9.00	-2.25
-10.25	-9.50	sph	-9.00	-0.75	-9.00	-0.75	-9.00	-0.75	-9.00	-1.25	-9.00	-1.25	-9.00	-1.25	-9.00	-1.75	-9.00	-1.75	-9.00	-2.25	-9.00	-2.25	-9.00	-2.25
-10.00	-9.00	sph	-9.00	-0.75	-9.00	-0.75	-9.00	-0.75	-9.00	-1.25	-9.00	-1.25	-9.00	-1.25	-9.00	-1.75	-9.00	-1.75	-9.00	-2.25	-9.00	-2.25	-9.00	-2.25
-9.75	-9.00	sph	-8.50	-0.75	-8.50	-0.75	-8.50	-1.25	-8.50	-1.25	-8.50	-1.25	-8.50	-1.25	-8.50	-1.75	-8.50	-1.75	-8.50	-2.25	-8.50	-2.25	-9.00	-2.25
-9.50	-8.50	sph	-8.50	-0.75	-8.50	-0.75	-8.50	-1.25	-8.50	-1.25	-8.50	-1.25	-8.50	-1.25	-8.50	-1.75	-8.50	-1.75	-8.50	-2.25	-8.50	-2.25	-8.50	-2.25
-9.25	-8.50	sph	-8.50	-0.75	-8.50	-0.75	-8.50	-1.25	-8.50	-1.25	-8.50	-1.25	-8.50	-1.25	-8.50	-1.75	-8.50	-1.75	-8.50	-2.25	-8.50	-2.25	8.50	-2.25
-9.00	-8.00	sph	-8.00	-0.75	-8.00	-0.75	-8.00	-1.25	-8.00	-1.25	-8.00	-1.25	-8.00	-1.25	-8.00	-1.75	-8.00	-1.75	-8.00	-2.25	-8.00	-2.25	-8.00	-2.25
-8.75	-8.00	sph	-8.00	-0.75	-8.00	-0.75	-8.00	-1.25	-8.00	-1.25	-8.00	-1.25	-8.00	-1.25	-8.00	-1.75	-8.00	-1.75	-8.00	-2.25	-8.00	-2.25	-8.00	-2.25
-8.50	-8.00	sph	-7.50	-0.75	-7.50	-0.75	-7.50	-1.25	-7.50	-1.25	-7.50	-1.25	-7.50	-1.25	-7.50	-1.75	-7.50	-1.75	-7.50	-2.25	-7.50	-2.25	-8.00	-2.25
-8.25	-7.75	sph	-7.50	-0.75	-7.50	-0.75	-7.50	-1.25	-7.50	-1.25	-7.50	-1.25	-7.50	-1.25	-7.50	-1.75	-7.50	-1.75	-7.50	-2.25	-7.50	-2.25	-7.50	-2.25
-8.00	-7.75	sph	-7.50	-0.75	-7.50	-0.75	-7.50	-1.25	-7.50	-1.25	-7.50	-1.25	-7.50	-1.25	-7.50	-1.75	-7.50	-1.75	-7.50	-2.25	-7.50	-2.25	-7.50	-2.25
-7.75	-7.50	sph	-7.00	-0.75	-7.00	-0.75	-7.00	-1.25	-7.00	-1.25	-7.00	-1.25	-7.00	-1.25	-7.00	-1.75	-7.00	-1.75	-7.00	-2.25	-7.00	-2.25	-7.00	-2.25
-7.50	-7.25	sph	-7.00	-0.75	-7.00	-0.75	-7.00	-1.25	-7.00	-1.25	-7.00	-1.25	-7.00	-1.25	-7.00	-1.75	-7.00	-1.75	-7.00	-2.25	-7.00	-2.25	-7.00	-2.25
-7.25	-7.00	sph	-6.50	-0.75	-6.50	-0.75	-6.50	-1.25	-6.50	-1.25	-6.50	-1.25	-6.50	-1.25	-6.50	-1.75	-6.50	-1.75	-6.50	-2.25	-6.50	-2.25	-6.50	-2.25
-7.00	-6.75	sph	-6.50	-0.75	-6.50	-0.75	-6.50	-1.25	-6.50	-1.25	-6.50	-1.25	-6.50	-1.25	-6.50	-1.75	-6.50	-1.75	-6.50	-2.25	-6.50	-2.25	-6.50	-2.25
-6.75	-6.50	sph	-6.00	-0.75	-6.00	-0.75	-6.00	-1.25	-6.00	-1.25	-6.00	-1.25	-6.00	-1.25	-6.00	-1.75	-6.00	-1.75	-6.00	-2.25	-6.00	-2.25	-6.00	-2.25
-6.50	-6.00	-0.75	-6.00	-0.75	-6.00	-0.75	-6.00	-1.25	-6.00	-1.25	-6.00	-1.25	-6.00	-1.25	-6.00	-1.75	-6.00	-1.75	-6.00	-2.25	-6.00	-2.25	-6.00	-2.25
-6.25	-5.75	-0.75	-5.75	-0.75	-5.75	-0.75	-5.75	-1.25	-5.75	-1.25	-5.75	-1.25	-5.75	-1.25	-5.75	-1.75	-5.75	-1.75	-5.75	-2.25	-5.75	-2.25	-5.50	-2.25
-6.00	-5.50	-0.75	-5.50	-0.75	-5.50	-0.75	-5.50	-1.25	-5.50	-1.25	-5.50	-1.25	-5.50	-1.25	-5.50	-1.75	-5.50	-1.75	-5.50	-2.25	-5.50	-2.25	-5.50	-2.25
-5.75	-5.50	-0.75	-5.50	-0.75	-5.50	-0.75	-5.50	-1.25	-5.50	-1.25	-5.50	-1.25	-5.50	-1.25	-5.50	-1.75	-5.50	-1.75	-5.50	-2.25	-5.25	-2.25	-5.50	-2.25
-5.50	-5.25	-0.75	-5.25	-0.75	-5.25	-0.75	-5.25	-1.25	-5.25	-1.25	-5.25	-1.25	-5.25	-1.25	-5.25	-1.75	-5.25	-1.75	-5.25	-2.25	-5.00	-2.25	-5.00	-2.25
-5.25	-5.00	-0.75	-5.00	-0.75	-5.00	-0.75	-5.00	-1.25	-5.00	-1.25	-5.00	-1.25	-5.00	-1.25	-5.00	-1.75	-5.00	-1.75	-5.00	-2.25	-5.00	-2.25	-5.00	-2.25
-5.00	-4.75	-0.75	-4.75	-0.75	-4.75	-0.75	-4.75	-1.25	-4.75	-1.25	-4.75	-1.25	-4.75	-1.25	-4.75	-1.75	-5.00	-1.75	-4.75	-2.25	-4.75	-2.25	-4.50	-2.25
-4.75	-4.50	-0.75	-4.50	-0.75	-4.50	-0.75	-4.50	-1.25	-4.50	-1.25	-4.50	-1.25	-4.50	-1.25	-4.75	-1.75	-4.50	-2.25	-4.50	-2.25	-4.50	-2.25		
-4.50	-4.25	-0.75	-4.25	-0.75	-4.50	-0.75	-4.25	-1.25	-4.25	-1.25	-4.25	-1.25	-4.25	-1.25	-4.25	-1.75	-4.50	-1.75	-4.25	-2.25	-4.25	-2.25		
-4.25	-4.00	-0.75	-4.00	-0.75	-4.25	-0.75	-4.00	-1.25	-4.00	-1.25	-4.00	-1.25	-4.00	-1.25	-4.00	-1.75	-4.25	-1.75	-4.00	-2.25	-4.25	-2.25		
-4.00	-4.00	-0.75	-4.00	-0.75	-4.00	-0.75	-4.00	-1.25	-4.00	-1.25	-4.00	-1.25	-4.00	-1.25	-4.00	-1.75	-4.25	-1.75	-4.00	-2.25	-4.00	-2.25		
-3.75	-3.75	-0.75	-3.75	-0.75	-3.75	-0.75	-3.75	-1.25	-3.75	-1.25	-3.75	-1.25	-3.75	-1.25	-3.75	-1.75	-3.75	-1.75	-3.75	-2.25	-4.00	-2.25		
-3.50	-3.50	-0.75	-3.50	-0.75	-3.50	-0.75	-3.50	-1.25	-3.50	-1.25	-3.50	-1.25	-3.50	-1.25	-3.50	-1.75	-3.50	-1.75	-3.50	-2.25	-3.75	-2.25		
-3.25	-3.25	-0.75	-3.25	-0.75	-3.25	-0.75	-3.25	-1.25	-3.25	-1.25	-3.25	-1.25	-3.25	-1.25	-3.25	-1.75	-3.25	-1.75	-3.25	-2.25	-3.50	-2.25		
-3.00	-3.00	-0.75	-3.00	-0.75	-3.00	-0.75	-3.00	-1.25	-3.00	-1.25	-3.00	-1.25	-3.00	-1.25	-3.00	-1.75	-3.00	-1.75	-3.00	-2.25	-3.25	-2.25		
-2.75	-2.75	-0.75	-2.75	-0.75	-2.75	-1.25	-2.75	-1.25	-2.75	-1.25	-2.75	-1.25	-2.75	-1.25	-2.75	-1.75	-2.75	-1.75	-2.75	-2.25	-3.00	-2.25		
-2.50	-2.50	-0.75	-2.50	-0.75	-2.50	-1.25	-2.50	-1.25	-2.50	-1.25	-2.50	-1.25	-2.50	-1.25	-2.50	-1.75	-2.75	-1.75	-2.50	-2.25	-2.75	-2.25		
-2.25	-2.25	-0.75	-2.25	-0.75	-2.25	-1.25	-2.25	-1.25	-2.25	-1.25	-2.25	-1.25	-2.25	-1.25	-2.25	-1.75	-2.50	-1.75	-2.25	-2.25	-2.50	-2.25		
-2.00	-2.00	-0.75	-2.00	-0.75	-2.00	-1.25	-2.00	-1.25	-2.00	-1.25	-2.00	-1.25	-2.00	-1.25	-2.00	-1.75	-2.25	-1.75	-2.00	-2.25	-2.25	-2.25		
-1.75	-1.75	-0.75	-1.75	-0.75	-1.75	-1.25	-1.75	-1.25	-1.75	-1.25	-1.75	-1.25	-1.75	-1.25	-2.00	-1.75	-2.00	-1.75	-1.75	-2.25	-2.00	-2.25		
-1.50	-1.50	-0.75	-1.50	-0.75	-1.50	-1.25	-1.50	-1.25	-1.50	-1.25	-1.50	-1.25	-1.50	-1.25	-1.75	-1.75	-1.50	-2.25	-1.50	-2.25	-1.75	-2.25		
-1.25	-1.25	-0.75	-1.25	-0.75	-1.25	-1.25	-1.25	-1.25	-1.25	-1.25	-1.25	-1.25	-1.75	-1.25	-1.50	-1.75	-1.25	-2.25	-1.50	-2.25	-1.50	-2.25		
-1.00	-1.00	-0.75	-1.00	-0.75	-1.00	-1.25	-1.00	-1.25	-1.00	-1.25	-1.00	-1.25	-1.00	-1.25	-1.25	-1.75	-1.00	-2.25	-1.25	-2.25	-1.25	-2.25		
-0.75	-0.75	-0.75	-0.75	-0.75	-0.75	-1.25	-0.75	-1.25	-0.75	-1.25	-0.75	-1.25	-0.75	-1.25	-1.00	-1.75	-0.75	-2.25	-1.00	-2.25	-1.00	-2.25		
-0.50	-0.50	-0.75	-0.50	-0.75	-0.50	-1.25	-0.50	-1.25	-0.50	-1.25	-0.50	-1.25	-0.50	-1.25	-0.75	-1.75	-0.50	-2.25	-0.75	-2.25	-0.75	-2.25		
-0.25	-0.25	-0.75	-0.25	-0.75	-0.25	-1.25	-0.25	-1.25	-0.25	-1.25	-0.25	-1.25	-0.25	-1.25	-0.50	-1.75	-0.25	-2.25	-0.50	-2.25	-0.50	-2.25		
+0.00	+0.00	-0.75	+0.00	-0.75	+0.00	-1.25	+0.00	-1.25	+0.00	-1.25	+0.00	-1.25	+0.00	-1.25	-0.25	-1.75	+0.00	-2.25	-0.25	-2.25	-0.25	-2.25		
+0.25	+0.25	-0.75	+0.25	-0.75	+0.25	-1.25	+0.25	-1.25	+0.25	-1.25	+0.25	-1.25	+0.25	-1.75	0.00	-1.75	+0.25	-2.25	+0.00	-2.25	+0.00	-2.25		
+0.50	+0.50	-0.75	+0.50	-0.75	+0.50	-1.25	+0.50	-1.25	+0.50	-1.25	+0.50	-1.75	+0.50	-1.75	+0.25	-1.75	+0.50	-2.25	+0.25	-2.25	+0.25	-2.25		
+0.75	+0.75	-0.75	+0.75	-0.75	+0.75	-1.25	+0.75	-1.25	+0.75	-1.25	+0.75	-1.75	+0.75	-1.75	+0.50	-1.75	+0.75	-2.25	+0.50	-2.25	+0.50	-2.25		
+1.00	+1.00	-0.75	+1.00	-0.75	+1.00	-1.25	+1.00	-1.25	+1.00	-1.25	+1.00	-1.75	+1.00	-1.75	+0.75	-1.75	+1.00	-2.25	+0.75	-2.25	+0.75	-2.25		
+1.25	+1.25	-0.75	+1.25	-0.75	+1.25	-1.25	+1.25	-1.25	+1.25	-1.25	+1.25	-1.75	+1.25	-1.75	+1.00	-1.75	+1.25	-2.25	+1.00	-2.25	+1.00	-2.25		
+1.50	+1.50	-0.75	+1.50	-0.75	+1.50	-1.25	+1.50	-1.25	+1.50	-1.25	+1.50	-1.75	+1.50	-1.75	+1.25	-1.75	+1.50	-2.25	+1.25	-2.25	+1.25	-2.25		
+1.75	+1.75	-0.75	+1.75	-0.75	+1.75	-1.25	+1.75	-1.25	+1.75	-1.25	+1.75	-1.75	+1.75	-1.75	+1.50	-1.75	+1.75	-2.25	+1.50	-2.25	+1.50	-2.25		

PARAMETRII LENTILĂ DE CONTACT DE PROBĂ



SCANEAZĂ
AICI



Întrebați reprezentantul Alcon sau vizitați hubul Alcon Fitting.

*pacienți astigmati cu cilindru de -2.50 sau mai mic
Pentru informații suplimentare referitoare la precauții, contraindicații și efecte adverse consultați instrucțiunile de utilizare.
©2022 Alcon RO-T3A-2300001

Alcon